

Control at your fingertips

Choosing a Riser Recliner


Independent Dual-motor

The independent dual-motor Riser Recliners (Medina) enable you to independently control the backrest and footrest using an easy 5-button power on handset – you can choose to raise the footrest and recline the backrest, or stay in a sitting position with only the footrest raised.



Multi-functional Dual-motor

The multi-functional dual-motor Ambassador Riser Recliners enable you to independently control the backrest and footrest using an easy to operate 6-button handset – you can choose to raise the footrest and recline the backrest, or stay in a sitting position with only the footrest raised. The smooth lift and tilt action also makes getting up from a seated position so much easier.

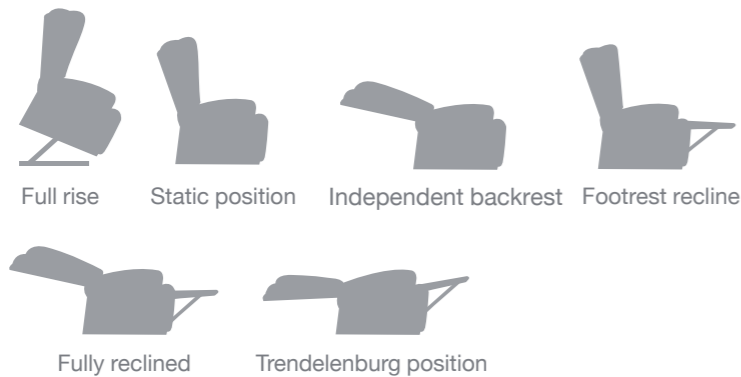

5-button handset
Independent
Dual-motor
Riser Recliners

6-button handset
Multi-functional
Dual-motor
Riser Recliners


Single-motor

The single-motor Riser Recliner (Lilburn) allows you to control the back rest and footrest at the same time using a simple to operate 3-button power on handset – the footrest automatically rises to support your legs when the back is reclined.




Single-motor Tilt-in-Space

Even when fully reclined, the Tilt-in-Space (Ellen) design retains the same angle between the seat and backrest, resulting in a higher raised feet position.

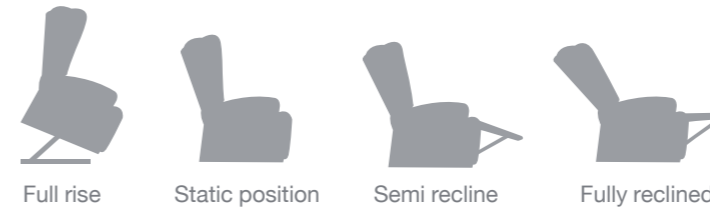

Single-motor Wallhugger

The clever design of the Wallhugger Riser Recliner (Walden) means you can position it close to a wall, even when fully reclined – as the back reclines, the seat base moves forward automatically and the foot rest rises.

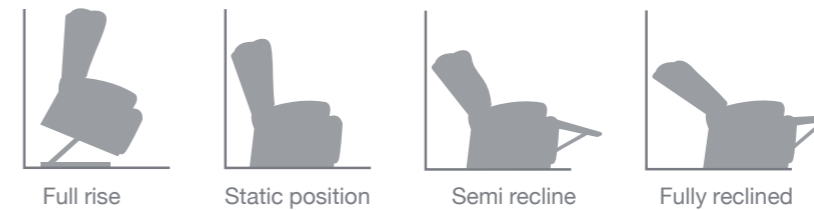

3-button handset
Single-motor
Riser Recliners